

● Analysis Report

Bruker IVDr Quantification in Plasma/Serum B.I.Quant-PS™

Sample ID: P_770_expno10.100000.11r

Measuring Date: 27-Mar-2023 18:56:41

Reporting Date: 27-Mar-2023 21:14:31, 7 page(s), Version 2.0.0

Quantification Method Version: Quant-PS 2.0.0

Disclaimer

RESEARCH USE ONLY: This is no clinical diagnostic analysis report. Must not be used for clinical (medical or IVD) diagnosis or for patient management! Additional concentration range information (95% range) provided numerically or graphically in this report must not be used for clinical diagnostic interpretation.

Application of B.I.QuantPS 2.0 requires use of Bruker's B.I.Methods SOP for plasma and serum.

Summary

The spectroscopic fingerprint of the sample is consistent with an EDTA-plasma profile.

The following metabolites were found with concentrations outside the 95% range of Bruker Quant-PS 2.0.0 plasma/serum metabolite concentration data base:

Amino acids and derivatives: Alanine (0.17 mmol/L), Histidine (0.04 mmol/L),

Carboxylic acids: 2-Hydroxybutyric acid (0.43 mmol/L), Formic acid (0.06 mmol/L), Lactic acid (1.61 mmol/L),

Keto acids and derivatives: 2-Oxoglutaric acid (0.06 mmol/L), Acetoacetic acid (0.06 mmol/L), Pyruvic acid (0.16 mmol/L),

Sugars and derivatives: Glucose (6.32 mmol/L),


Technical additives: K-EDTA (7.59 mmol/L).

Further detailed information is provided on the following pages.

Contents


1	Alcohols and derivatives	3
2	Amines and derivatives	3
3	Amino acids and derivatives	3
4	Carboxylic acids	4
5	Essential nutrient	4
6	Keto acids and derivatives	4
7	Sugars and derivatives	4
8	Sulfones	5
9	Technical additives	5
10	Explanations	6
10.1	How to read the result	6
10.1.1	Result parameters	6
10.1.2	Different fit situations	7
10.2	Detailed definitions	7

1 Alcohols and derivatives

Compound	Conc. mmol/L	LOD mmol/L	r mmol/L	ρ %	Δ mmol/L	95% Range mmol/L	Graphics (*)
Ethanol	< 0.10	0.10	0.000	0 ○	2.683	≤ 0.82	


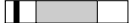



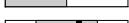













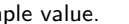

(*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

2 Amines and derivatives

Compound	Conc. mmol/L	LOD mmol/L	r mmol/L	ρ %	Δ mmol/L	95% Range mmol/L	Graphics (*)
Trimethylamine-N-oxide	< 0.08	0.08	0.041	94 ○	0.075	≤ 0.08	







(*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

3 Amino acids and derivatives

Compound	Conc. mmol/L	LOD mmol/L	r mmol/L	ρ %	Δ mmol/L	95% Range mmol/L	Graphics (*)
2-Aminobutyric acid	< 0.05	0.05	0.028	0 ○	0.142	≤ 0.10	
Alanine	0.17	0.02	0.175	99 ●	0.009	0.29 - 0.64	
Asparagine	< 0.05	0.05	0.000	0 ○	3.585	≤ 0.08	
Creatine	0.02	0.01	0.021	99 ●	0.003	≤ 0.07	
Creatinine	0.06	0.01	0.064	99 ●	0.003	0.06 - 0.14	
Glutamic acid	< 0.05	0.05	0.000	0 ○	1.117	≤ 0.24	
Glutamine	0.67	0.02	0.668	99 ●	0.026	0.30 - 0.83	
Glycine	0.25	0.01	0.249	100 ●	0.006	0.17 - 0.44	
Histidine	0.04	0.02	0.040	99 ●	0.002	0.07 - 0.16	
Isoleucine	0.04	0.03	0.044	99 ●	0.004	0.03 - 0.11	
Leucine	0.08	0.01	0.077	98 ●	0.008	0.07 - 0.20	
Lysine	0.22	0.04	0.221	73 ○	0.061	≤ 0.29	
Methionine	< 0.05	0.05	0.011	91 ○	0.002	0.05 - 0.13	
N,N-Dimethylglycine	< 0.01	0.01	0.004	68 ○	0.001	≤ 0.01	
Ornithine	0.04	0.02	0.038	96 ●	0.017	≤ 0.16	
Phenylalanine	0.07	0.03	0.069	98 ●	0.004	≤ 0.07	
Proline	0.31	0.05	0.310	43 ○	0.147	≤ 0.59	
Sarcosine	< 0.01	0.01	0.008	96 ●	0.001	≤ 0.01	
Threonine	0.08	0.04	0.077	64 ○	0.054	≤ 0.24	
Tyrosine	0.07	0.03	0.069	97 ●	0.005	≤ 0.08	
Valine	0.24	0.03	0.236	100 ●	0.004	0.15 - 0.35	


(*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

4 Carboxylic acids

Compound	Conc. mmol/L	LOD mmol/L	r mmol/L	ρ %	Δ mmol/L	95% Range mmol/L	Graphics (*)
2-Hydroxybutyric acid	0.43	0.15	0.428	58 ○	0.121	≤ 0.17	
Acetic acid	0.04	0.01	0.036	99 ●	0.001	≤ 0.06	
Citric acid	0.11	0.03	0.110	95 ●	0.020	≤ 0.21	
Formic acid	0.06	0.02	0.060	100 ●	0.001	≤ 0.03	
Lactic acid	1.6	0.03	1.614	100 ●	0.046	2.23 - 7.14	
Succinic acid	< 0.01	0.01	0.005	99 ●	0.000	≤ 0.01	






(*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

5 Essential nutrient

Compound	Conc. mmol/L	LOD mmol/L	r mmol/L	ρ %	Δ mmol/L	95% Range mmol/L	Graphics (*)
Choline	< 0.05	0.05	0.000	0 ○	0.111	≤ 0.06	

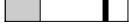


(*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

6 Keto acids and derivatives

Compound	Conc. mmol/L	LOD mmol/L	r mmol/L	ρ %	Δ mmol/L	95% Range mmol/L	Graphics (*)
2-Oxoglutaric acid	0.06	0.02	0.058	78 ○	0.018	≤ 0.02	
3-Hydroxybutyric acid	0.20	0.02	0.199	97 ●	0.030	≤ 0.26	
Acetoacetic acid	0.06	0.01	0.063	99 ●	0.002	≤ 0.02	
Acetone	0.05	0.01	0.050	99 ●	0.001	≤ 0.06	
Pyruvic acid	0.16	0.03	0.156	100 ●	0.001	≤ 0.07	


(*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

7 Sugars and derivatives

Compound	Conc. mmol/L	LOD mmol/L	r mmol/L	ρ %	Δ mmol/L	95% Range mmol/L	Graphics (*)
D-Galactose	< 0.11	0.11	0.000	0 ○	2.153	≤ 0.11	
Glucose	6.3	0.54	6.320	100 ●	0.018	1.73 - 6.08	
Glycerol	< 0.08	0.08	0.042	0 ○	0.271	≤ 0.44	



(*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

8 Sulfones

Compound	Conc. mmol/L	LOD mmol/L	r mmol/L	ρ %	Δ mmol/L	95% Range mmol/L	Graphics (*)
Dimethylsulfone	< 0.01	0.01	0.009	99 ●	0.000	≤ 0.02	

(*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

9 Technical additives

Compound	Conc. mmol/L	LOD mmol/L	r mmol/L	ρ %	Δ mmol/L	95% Range mmol/L	Graphics (*)
Ca-EDTA	2.1	0.50	2.148	100 ●	0.009	1.90 - 2.60	
K-EDTA	7.6	0.50	7.591	100 ●	0.032	4.50 - 7.50	

(*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

10 Explanations

This section contains the definition of the parameters used above. In the section 10.1 a short manual, how to interpret the results, is presented. The section 10.2 contains the exact definitions of the parameters r , ρ and Δ .

10.1 How to read the result



Figure 1: Examples of fitting.

In the figure 1(a), the black line, the blue line and the yellow line represent the original spectrum, the calculated signal fit and its baseline, respectively.

The blue area relates to the metabolite concentration to be determined and the red area represents a residue.

In case of the signal overlap a different approach is used: two or more overlapping signals are being fitted simultaneously. The most iconic example of such signals are the ones generated by CH_3 groups of Creatinine and Creatine. In such a case, the blue line and the grey area relate the sum of all fitted signals. The blue area corresponds to the concentration of the metabolite of interest (cf. figure 1(b)).

10.1.1 Result parameters

- Conc.** is the final result concentration of the metabolite,
- LOD** is the *limit of detection* of the given metabolite,
- r is the *raw concentration* i.e. the concentration equivalent of the resulting signal fit prior to comparing to **LOD** (relates to the blue area, cf. α),
- ρ is the correlation of lineshape metabolite signal with calculated fit characterizing the match between metabolite signal and fit (cf. β),
- Δ is the concentration equivalent of the difference between metabolite signal and calculated fit (residue corresponding to the red area, cf. γ).

10.1.2 Different fit situations

Now we will describe the main fit cases.

- i) In an ideal situation, where the fit corresponds fully to a metabolite signal well above **LOD**:
 - the raw concentration is similar to the final result concentration,
 - the correlation is $\geq 95\%$ (indicated by ● displayed next to the value, otherwise the mark ○ is being used),
 - the residue Δ is close to zero mmol/L.
- ii) Similar to situation described in i), but raw concentration below **LOD**. Generally, only an upper limit (e.g. $< \text{LOD}$) can be provided. Especially, if the difference between raw concentration **r** and the final concentration **Conc.** is small, use the graphical figure displayed when pointing the cursor on the metabolite name for further visual inspection and validation. If a metabolite signal can be clearly discriminated from the rest of the spectrum and the calculated fit represents the respective signal well, the raw concentration may be used as approximative concentration estimate.
- iii) Low correlation combined with large residual Δ . Such situation may arise in case of significant signal overlap, e.g. if a doublet signal of a metabolite to be quantified is overlaid with a large singlet. Use the graphical figure displayed when pointing the cursor on the metabolite name for further visual inspection and validation. If the non-overlaid part of the signal is well fitted, the final calculated concentration may still be used with confidence depending on the degree and nature of signal overlap.
- iv) Combination of ii) and iii). Use the graphical figure displayed when pointing the cursor on the metabolite name for further visual inspection and validation.

10.2 Detailed definitions

Let s , f and b denote the functions describing the *raw spectra*, *fitted curve* and *(fitted) baseline* respectively. These functions are chosen such that $s \approx f + b$. Moreover, let I be a relevant PPM interval and P_N be the proton number for given metabolite/signal.

α) **r** (*raw concentration*) is defined as

$$\mathbf{r} = \frac{1}{P_N} \int_{\mathbb{R}} f(\xi) \, d\xi.$$

β) ρ is the *correlation* of the functions s and $f + b$, i.e.

$$\rho = \max(0, \text{corr}(\bar{s}, \overline{f + b})),$$

where \bar{s} , $\overline{f + b}$ are numerical representations of the functions s and $f + b$ on sufficiently fine mesh of the interval I .

γ) Δ is the the area between the raw signal s and the fitted data $f + b$ on the interval I expressed in the term of the concentration, i.e.

$$\Delta = \frac{1}{P_N} \int_I |s(\xi) - f(\xi) - b(\xi)| \, d\xi.$$